
Support wire guides through challenging anatomy with the catheter that has a crafted tip.

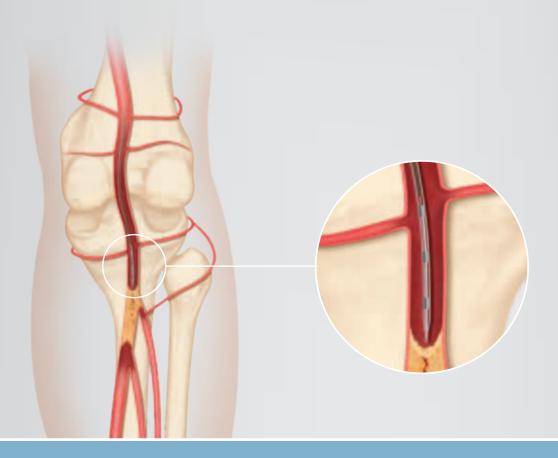

CrossCath[®]

SUPPORT CATHETER

Crafted Tip

CrossCath® features a tip design that transitions smoothly to the wire and is engineered to help cross lesions and minimize the risk of vessel trauma.

Excellent Pushability

The tapered shaft enhances pushability through tortuous anatomy and challenging lesions.

Low Profile

A distal end as small as 1.9 Fr allows access to smaller vasculature.





Enlarged view of the distal tip.

CrossCath^{*}

CrossCath Support Catheters are single-lumen intravascular catheters designed to support a wire guide during access of vasculature, allow for the exchange of wire guides, and provide a conduit for the delivery of saline solutions or diagnostic contrast agents. CrossCath also has these features:

- Three highly radiopaque platinum-iridium markers that aid in estimating lengths within the vasculature.
- A hydrophilic coating on the distal end for smooth advancement.

| Order Number | Reference Part Number | Proximal/Distal Fr | Accepts Wire Guide Diameter inch | Length cm | Tip Configuration |
|--|--------------------------|-----------------------|---|--------------|----------------------|
| With Hydrophilic Coating on Distal 40 cm | | | | | |
| G18495 | CXC3.0-1.9-14-65-P-NS-0 | 3.0/1.9 | .014 | 65 | straight |
| G51572 | CXC3.0-1.9-14-90-P-NS-0 | 3.0/1.9 | .014 | 90 | straight |
| G51573 | CXC3.0-1.9-14-135-P-NS-0 | 3.0/1.9 | .014 | 135 | straight |
| G51574 | CXC3.0-1.9-14-150-P-NS-0 | 3.0/1.9 | .014 | 150 | straight |
| G18497 | CXC3.4-2.2-18-65-P-NS-0 | 3.4/2.2 | .018 | 65 | straight |
| G51575 | CXC3.4-2.2-18-90-P-NS-0 | 3.4/2.2 | .018 | 90 | straight |
| G51583 | CXC3.4-2.2-18-135-P-NS-0 | 3.4/2.2 | .018 | 135 | straight |
| G51577 | CXC3.4-2.2-18-150-P-NS-0 | 3.4/2.2 | .018 | 150 | straight |
| G18499 | CXC4.8-3.7-35-65-P-NS-0 | 4.8/3.7 | .035 | 65 | straight |
| With Hydrophilic Coating on Distal 60 cm | | | | | |
| G51580 | CXC4.8-3.7-35-90-P-NS-0 | 4.8/3.7 | .035 | 90 | straight |
| G51581 | CXC4.8-3.7-35-135-P-NS-0 | 4.8/3.7 | .035 | 135 | straight |
| G51582 | CXC4.8-3.7-35-150-P-NS-0 | 4.8/3.7 | .035 | 150 | straight |

Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Service for details.

Customer Service

EMEA: EDI - www.cookmedical.com/edi.do

Distributors: +353 61239240, ssc.distributors@cookmedical.com Austria: +43 179567121, oe.orders@cookmedical.com

Belgium: +32 27001633, be.orders@cookmedical.com Denmark: +45 38487607, da.orders@cookmedical.com Finland: +358 972519996, fi.orders@cookmedical.com France: +33 171230269, fr.orders@cookmedical.com Germany: +49 6950072804, de.orders@cookmedical.com Hungary: +36 17779199, hu.orders@cookmedical.com Ireland: +353 61239252, ie.orders@cookmedical.com Italy: +39 0269682853, it.orders@cookmedical.com

Netherlands: +31 202013367, nl.orders@cookmedical.com Norway: +47 23162968, no.orders@cookmedical.com Spain: +34 912702691, es.orders@cookmedical.com Sweden: +46 858769468, se.orders@cookmedical.com

Switzerland - French: +41 448009609, fr.orders@cookmedical.com Switzerland - Italian: +41 448009609, it.orders@cookmedical.com Switzerland - German: +41 448009609, de.orders@cookmedical.com United Kingdom: +44 2073654183, uk.orders@cookmedical.com

www.cookmedical.com

Americas: EDI - www.cookmedical.com/edi.do Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335 E-mail: orders@cookmedical.com

Australia:

Phone: +61 738411188, 1800777222, Fax: +61 738411288, 1800077283

E-mail: cau.custserv@cookmedical.com

