


## Perform two procedures with one device.



Sonohysterography & Biopsy

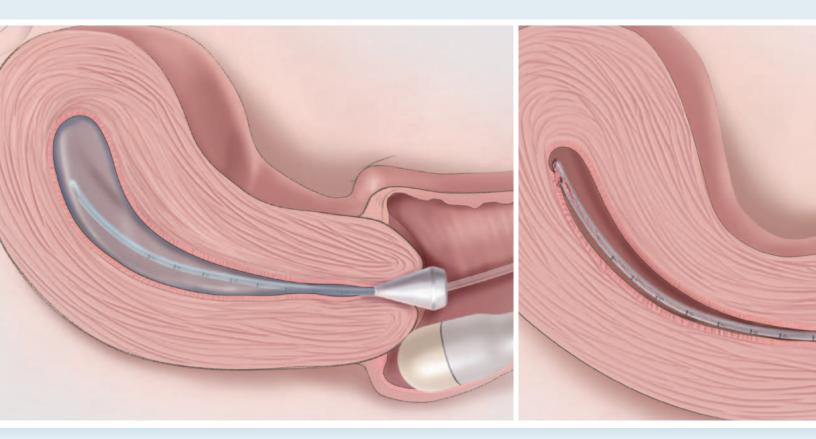

## Improve the visualization of the uterine cavity and obtain an endometrial biopsy.

The Goldstein SonoBiopsy Catheter allows physicians to perform both sonohysterography and endometrial biopsy in one procedure, sparing the patient the discomfort of two speculum placements. This specially designed catheter has a round, closed tip that makes it easier to access the uterus and has a wide sideport that facilitates a seamless transition between the instillation of saline solution and the aspiration of endometrial tissue.

Cook Medical-simplifying procedures for patients and clinicians.

## Sonohysterography



The Goldstein SonoBiopsy Catheter is conveniently marked in 1 cm graduations up to 10 cm to indicate proper positioning. The catheter also features an adjustable acorn positioner that fits snugly into the cervical os, minimizing leakage.

## Biopsy

Sideport effectively aspirates endometrial tissue.



10

COOK Goldstein SONOBIOPSY<sup>™</sup> CATHETER

Global Product Number	Order Number	Fr	Length cm	Unit of Sale
G46190	J-GSBX-072026	7.2	26	box of 10

Visit www.cookmedical.com/sb to view a step-by-step procedural video.



COOK MEDICAL INCORPORATED PO. Box 4195, Bloomington, IN 47402-4195 USA Phone: 812.339.2235, Toll Free: 800.457.4500, Toll Free Fax: 800.554.8335

COOK IRELAND LTD. O'Halloran Road, National Technology Park Limerick, IRELAND, Phone: +353 61 334440, Fax: +353 61 334441

www.cookmedical.com

COOK (CANADA) INC. 111 Sandiford Drive, Stouffville, Ontario, L4A 7X5 CANADA Phone: 905.640.7110, Toll Free: 800.668.0300, Fax: 905.640.6804

WILLIAM A. COOK AUSTRALIA PTY. LTD. 95 Brandl Street, Brisbane Technology Park, Eight Mile Plains Brisbane, QLD 4113 AUSTRALIA, Phone: +61 7 3841 1188

© COOK 2010

WH-BM-GSBX-EN-201003

AORTIC	CRITICAL CARE	ENDOSCOPY	INTERVENTIONAL RADIOLOGY	LEAD MANAGEMENT	PERIPHERAL	SURGERY	UROLOGY	WOMEN'S HEALTH	
	0/11/2		1010101001						