


Pack tapered Tornado[®] coils in selective or short vessels.

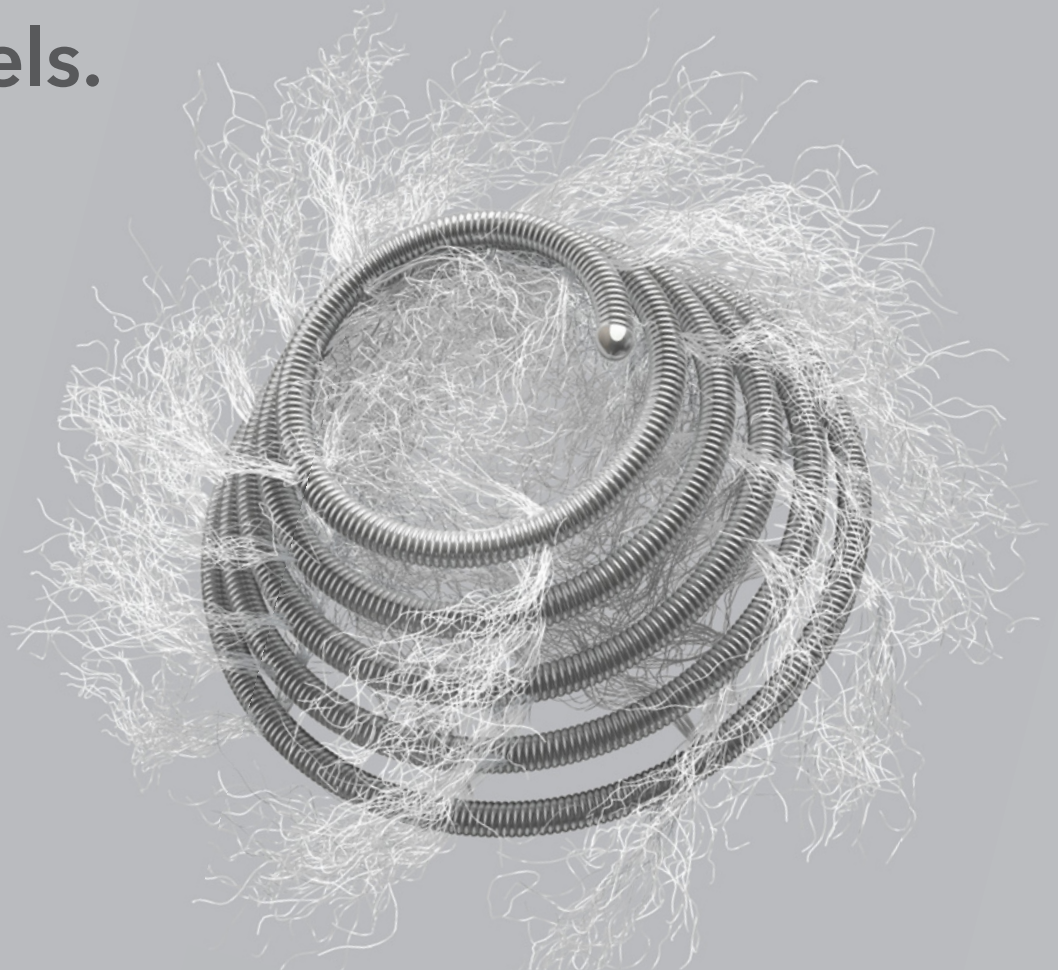

Tornado[®] EMBOLIZATION COIL

Ideal for tapered vessels

The Tornado is the choice for smaller vessels that taper down in width. The coil can also occlude short, straight vessels.

A large range of coil profiles


The Tornado has a large number of profiles that are designed for use in .018 and .035 inch coil systems.

Excellent packing coil

The soft platinum Tornado coil conforms to the vessel to create a tight, occluding mass.



www.cookmedical.com

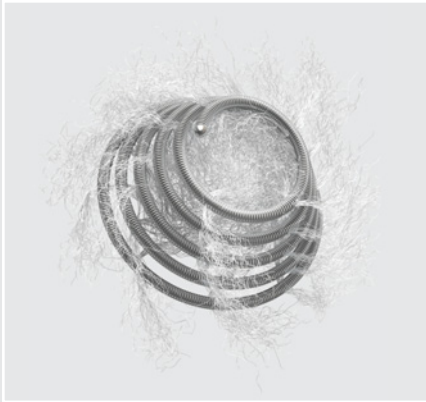

Tornado®

EMBOLIZATION COIL

The platinum coil with synthetic fibers is used for the embolization of selective vessel supply to arteriovenous malformations and other vascular lesions.

The tornado-shaped configuration maximizes coil exposure to a cross-section of the lumen for disruption of blood flow.

The soft platinum conforms to the vessel to create a tight occluding mass.

The coil is ideally suited for tapering vessel situations.

The coil is available in sizes that fit .018 and .035 inch catheters.

Tapering diameters range from 4-3 to 10-5 mm.

The coil exits the delivery catheter small end first.

Order Number	Reference Part Number	Recommended Catheter ID and End Hole Diameter inch	Extended Embolus Length cm	Coiled Embolus Tapering Diameter mm
.018 inch Microcoils				
G08261	MWCE-18S-3/2-TORNADO	.018	2	3 - 2
G08357	MWCE-18S-4/2-TORNADO	.018	4	4 - 2
G08356	MWCE-18S-5/2-TORNADO	.018	5	5 - 2
G08259	MWCE-18S-6/2-TORNADO	.018	7	6 - 2
G09218	MWCE-18S-7/3-TORNADO	.018	9	7 - 3
G09215	MWCE-18S-8/4-TORNADO	.018	9.5	8 - 4
G09217	MWCE-18S-10/4-TORNADO	.018	14.2	10 - 4
.035 inch Coils				
G10414	MWCE-35-4/3-TORNADO	.035	2.6	4 - 3
G10415	MWCE-35-5/3-TORNADO	.035	4.1	5 - 3
G10416	MWCE-35-6/3-TORNADO	.035	5.8	6 - 3
G10417	MWCE-35-7/3-TORNADO	.035	8	7 - 3
G10411	MWCE-35-8/4-TORNADO	.035	8	8 - 4
G10412	MWCE-35-8/5-TORNADO	.035	8.2	8 - 5
G10428	MWCE-35-10/4-TORNADO	.035	12.5	10 - 4
G10413	MWCE-35-10/5-TORNADO	.035	12.5	10 - 5

Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Service for details.

Customer Service

EMEA: EDI - www.cookmedical.com/edi.do
 Distributors: +353 61239240, ssc.distributors@cookmedical.com
 Austria: +43 179567121, oe.orders@cookmedical.com
 Belgium: +32 27001633, be.orders@cookmedical.com
 Denmark: +45 38487607, da.orders@cookmedical.com
 Finland: +358 972519996, fi.orders@cookmedical.com
 France: +33 171230269, fr.orders@cookmedical.com
 Germany: +49 6950072804, de.orders@cookmedical.com
 Hungary: +36 17779199, hu.orders@cookmedical.com
 Ireland: +353 61239252, ie.orders@cookmedical.com
 Italy: +39 0269682853, it.orders@cookmedical.com
 Netherlands: +31 202013367, nl.orders@cookmedical.com
 Norway: +47 23162968, no.orders@cookmedical.com
 Spain: +34 912702691, es.orders@cookmedical.com
 Sweden: +46 858769468, se.orders@cookmedical.com
 Switzerland - French: +41 448009609, fr.orders@cookmedical.com
 Switzerland - Italian: +41 448009609, it.orders@cookmedical.com
 Switzerland - German: +41 448009609, de.orders@cookmedical.com
 United Kingdom: +44 2073654183, uk.orders@cookmedical.com

www.cookmedical.com

Americas: EDI - www.cookmedical.com/edi.do
 Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335
 E-mail: orders@cookmedical.com

Australia:
 Phone: +61 738411188, 1800777222, Fax: +61 738411288, 1800077283
 E-mail: cau.custserv@cookmedical.com

