
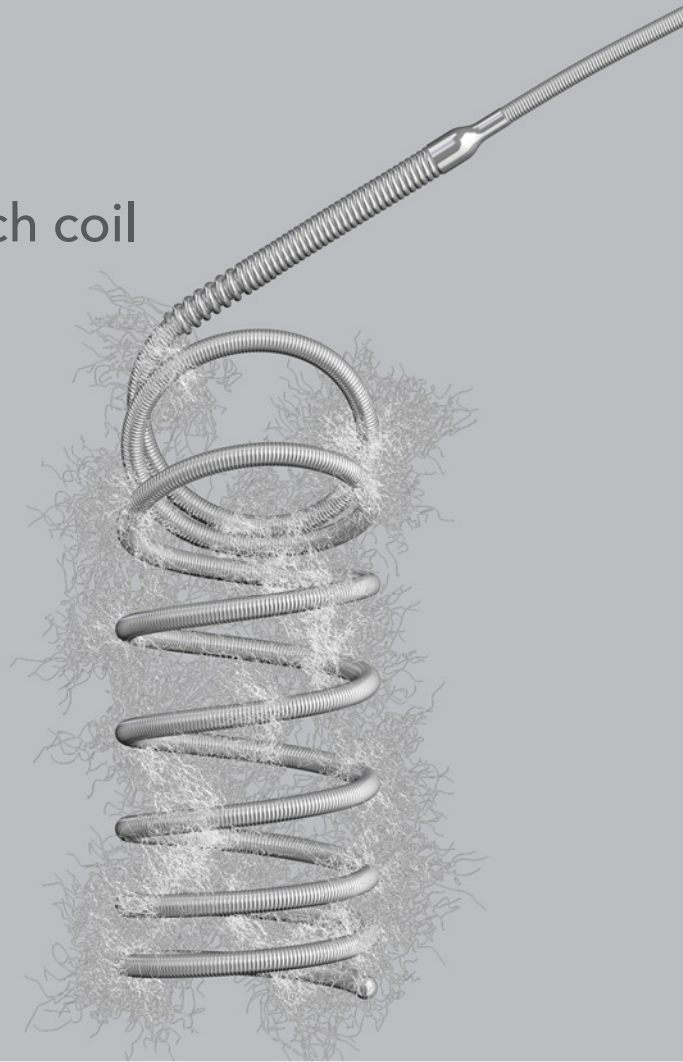


Choose simplicity.

Retracta™ is a detachable .035 inch coil that's fully retractable.

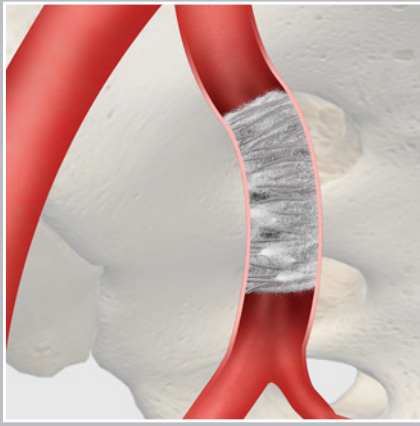

Retracta™

DETACHABLE EMBOLIZATION COIL

Use it right out of the package.

Retracta is easy to use and requires minimal preparation time. It's an economical choice that gives you another viable option for treating patients.

Use a detachable option based on platinum Nester® technology.


Retracta's soft platinum lets you pack the coils tightly. Use Retracta with the popular Nester and Tornado® coils to control the embolus size.

Detach with precision when you are ready.

Place the coil and detach it. Or if you need to change its placement, retract and reposition the coil before detaching it. Under fluoroscopy, the change in radiopacity at the detachment zone tells you where the coil is. Then you can inject contrast to confirm that the coil is correctly placed.



www.cookmedical.com

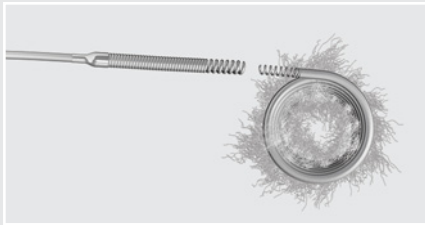

Retractable™

DETACHABLE EMBOLIZATION COIL

- Retractable Detachable Embolization Coils are intended for arterial and venous embolization in the peripheral vasculature.
- The coils can be advanced and fully repositioned before they are deployed.

Order Number	Reference Part Number	Recommended Catheter ID/End Hole Diameter inch	Extended Embolus Length cm	Coiled Embolus Diameter mm	Approximate Number of Loops
.035 Inch Detachable Coils					
G34188	MWCER-35-7-6	.035	7	6	3.7
G34189	MWCER-35-7-8	.035	7	8	2.8
G34190	MWCER-35-7-10	.035	7	10	2.2
G34191	MWCER-35-7-12	.035	7	12	1.9
G34192	MWCER-35-7-14	.035	7	14	1.6
G34193	MWCER-35-7-16	.035	7	16	1.4
G34194	MWCER-35-7-18	.035	7	18	1.2
G34195	MWCER-35-7-20	.035	7	20	1.1
G34197	MWCER-35-14-6	.035	14	6	7.4
G34198	MWCER-35-14-8	.035	14	8	5.6
G34199	MWCER-35-14-10	.035	14	10	4.5
G34200	MWCER-35-14-12	.035	14	12	3.7
G34201	MWCER-35-14-14	.035	14	14	3.2
G34202	MWCER-35-14-16	.035	14	16	2.8
G34203	MWCER-35-14-18	.035	14	18	2.5
G34204	MWCER-35-14-20	.035	14	20	2.2
Order Number	Reference Part Number	Accepts Catheter Fr	Accepts Instrumentation Diameter inch		
Tuohy-Borst Sidearm Adapter (optional accessory for contrast injection)					
G06299	PTBYC-RA	≤ 9.0	.014 thru .118		

Some products or part numbers may not be available in all markets. Contact your local Cook representative or Customer Service for details.

Customer Service

EMEA: EDI - www.cookmedical.com/edi.do
 Distributors: +353 61239240, ssc.distributors@cookmedical.com
 Austria: +43 179567121, oe.orders@cookmedical.com
 Belgium: +32 27001633, be.orders@cookmedical.com
 Denmark: +45 38487607, da.orders@cookmedical.com
 Finland: +358 972519996, fi.orders@cookmedical.com
 France: +33 171230269, fr.orders@cookmedical.com
 Germany: +49 6950072804, de.orders@cookmedical.com
 Hungary: +36 17779199, hu.orders@cookmedical.com
 Ireland: +353 61239252, ie.orders@cookmedical.com
 Italy: +39 0269682853, it.orders@cookmedical.com
 Netherlands: +31 202013367, nl.orders@cookmedical.com
 Norway: +47 23162968, no.orders@cookmedical.com
 Spain: +34 912702691, es.orders@cookmedical.com
 Sweden: +46 858769468, se.orders@cookmedical.com
 Switzerland - French: +41 448009609, fr.orders@cookmedical.com
 Switzerland - Italian: +41 448009609, it.orders@cookmedical.com
 Switzerland - German: +41 448009609, de.orders@cookmedical.com
 United Kingdom: +44 2073654183, uk.orders@cookmedical.com

www.cookmedical.com

Americas: EDI - www.cookmedical.com/edi.do
 Phone: +1 812.339.2235, 800.457.4500, Fax: 800.554.8335
 E-mail: orders@cookmedical.com

Australia:
 Phone: +61 738411188, 1800777222, Fax: +61 738411288, 1800077283
 E-mail: cau.custserv@cookmedical.com

